



SUPER COMPACT LABORATORY CENTRIFUGES RA 20 Vx



RA 20 Vx

KEY FEATURES OF RA TYPE CENTRIFUGES



EASE OF USE

- Compact and light equipment, easy to be moved.
- Sight glass.
- Front control panel with digital speed display.



ADAPTABILITY / PERFORMANCE

- Interchanging elements.
- Possibility to test various techniques (filtration / decantation) for liquid / solid separation.
- Ideal for industries such as Chemistry, Pharmacy, Food or Textile Industries.



FONCTIONNALITY

- Ideal for R&D departments.
- Compact and functional design.
- Lid clamping latches (for increased tightness).

BASIC UNIT FEATURES




MAIN FEATURES

TYPE	Ø Basket mm	Basket Height mm	Useful volume l	Useful load kg	Maximum speed rpm	Length mm	Width mm	Height mm	Weight kg
RA 20	200	100	1.2	1.5	3000	516	454	440	55
	120		0.6	0.75					
	280		2.75	3.45					
	120	50	0.3	0.37					



ALTERNATIVE CONSTRUCTIONS AND OPTIONAL EQUIPMENTS

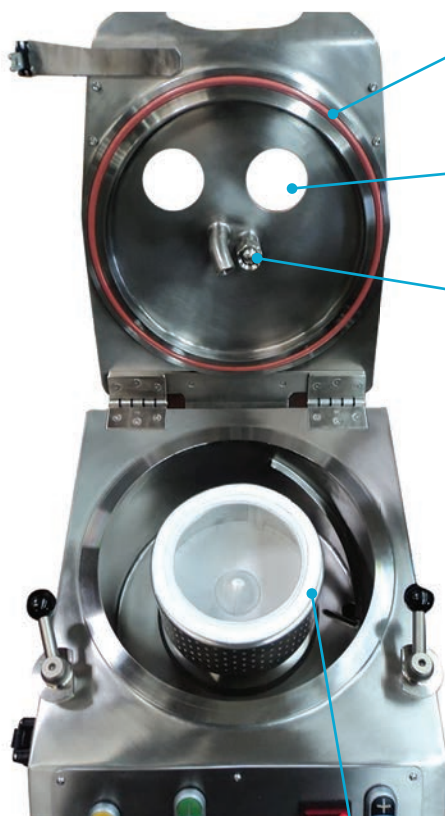


Lid gasket
in silicone rubber PTFE
encapsulated

Additional sight glass

Cake washing pipe with
orientable nozzle
(RA 20 Vx R).

Different diameter baskets
(120 and 280 mm) with useful volu-
mes from 0.3L to 2.75L (in accor-
dance with diameter and height).

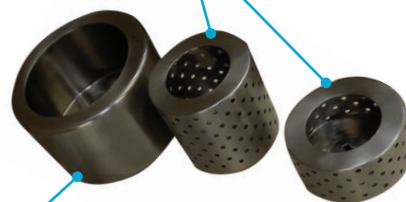


Decantation bowls
(purified liquid discharged by overflow)
with diameters of 120, 200 and 280 mm.

Possibility of designing the bowl with set(s) of
holes plugged with screws:

- In the top to modify capacity
- In the bottom plate to drain the liquid phase
once the basket has stopped.



Neutral gas feeding piping
allowing to work under
flowing nitrogen.

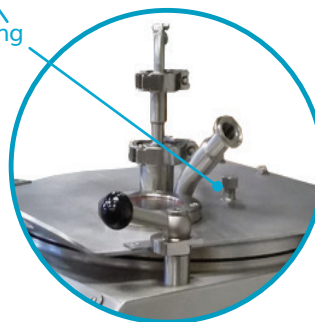

Removable filter bag
(different textile fabrics with various
features available) / Draining cloth.

OPTIONAL EQUIPMENT AND AUXILIARIES

Maximum speed increased to 3.500 or 4.000
rpm (in accordance with diameter, type of rotor
and features of products to be processed).

Alternative materials of construction: other stain-
less steel grades (AISI 904), Alloy C, coating of shell,
frame and lid...

**Various polishing levels, connection types and
voltage.**



Timer
Cycle time selection.

SPECIFIC VERSIONS

MANUAL TROLLEY / SKID

With wheels for the centrifuge only or with selected peripheral equipment.



On skid with additional basket.



With peripheral equipment (storage containers for feeding solution, cake washing liquid and filtrate, and general control cabinet mounted on manual trolley with wheels).

COATING MATERIALS

Manufactured from stainless steel HALAR coated and Alloy C, with gasket on lid, cake washing pipe, neutral gas feeding piping and vent piping on outlet pipe.



HEADQUARTERS

ROUSSELET ROBATEL
45 Avenue Rhin et Danube
Parc d'activité économique de Marenton
07104 ANNONAY - FRANCE
Tel.: +33 / (0) 4 75 69 22 11
Fax: +33 / (0) 4 75 67 69 80
E-mail: rousselet.sa@rousselet.com

UK

Rousselet UK Ltd
Parkside House, 17 East Parade
HARROGATE
NORTH YORKSHIRE HG1 5LF
Tel: + 44 (0)1 423 530 093
Fax: + 44 (0)1 423 530 120
E-mail: sales@rousselet-robatel.co.uk

CZECH REPUBLIC

Rousselet Robatel CZ
Počernická 96
CZ – 108 00 Praha 10
Tel: +420 296 411 402
Fax: +420 296 411 403
E-mail: info@rousselet.cz

USA

Robatel Inc.
703 West Housatonic Street
PITTSFIELD
MA 01201
Tel: + 1 413 499 4818
Fax: + 1 413 499 5648
E-mail: sales@robatel.com

WWW.ROUSSELET-ROBATEL.COM